

Omnitrope[®]
Somatropin

POWERED BY

SurePal[®]
Device

a Sandoz innovation




Information for patients:
HOW TO USE GUIDE

SurePal[®] 5


For use with Omnitrope® 5 mg/1.5 mL pre-assembled cartridges.
Remember that only Omnitrope® 5 mg/1.5 mL cartridges will fit in your SurePal® 5.
If the cartridges you have do not fit easily, contact your healthcare professional.

For full instructions, please refer to the SurePal® 5 User's Guide
supplied with your SurePal® pen.

Parts of your SurePal® 5

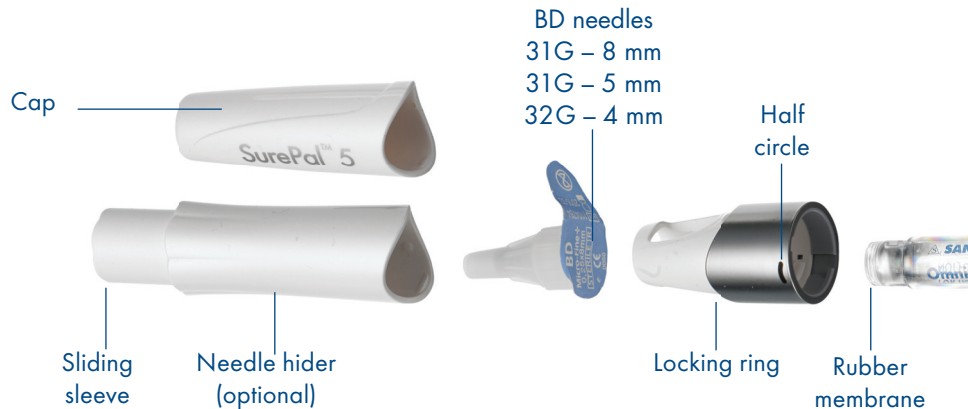

Cartridges must always be stored in the refrigerator, between 2 and 8 degrees Celsius. It is recommended that you take the pen and cartridge(s) out of the refrigerator for 30 minutes before your injection, and return to refrigerator after you inject.

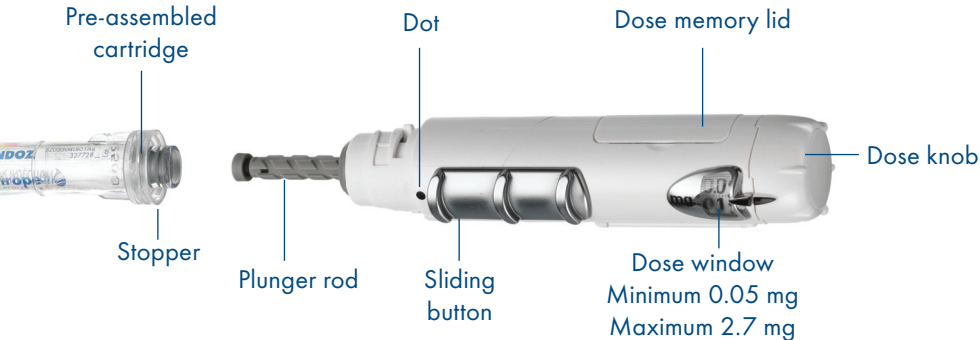


In addition to the items shown here, you will need a sharps box and disinfectant (alcohol swabs). Pre-assembled cartridges, alcohol swabs and needles do not come with your SurePal®.

1. How to prepare your SurePal[®] pen for your injections

Wash your hands. Begin with your SurePal[®] taken apart (as on the back page).

Set the dose memory (when the dose changes).


- Open the dose memory lid and set the dose by turning the dose setting wheel until the exact dose prescribed by your doctor is displayed and lines up with the lock symbol .
- In the below example the dose has been set to 2.4 mg.

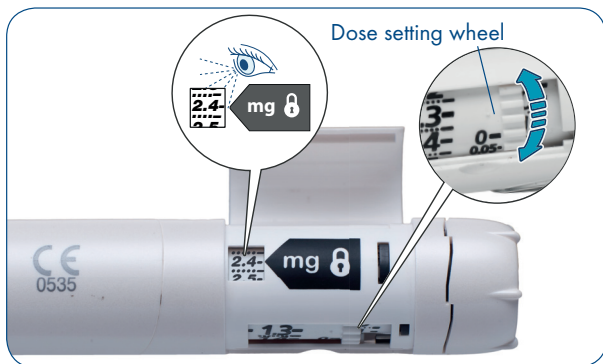

- Odd numbers are not shown but indicated by a dotted line “- -” (e.g. 2.5 mg is the dotted line between 2.4 mg and 2.6 mg.)
- Close the dose memory lid and ensure it clicks into place.

Attach the pre-assembled cartridge


- To insert a cartridge into the locking ring, hold the cartridge with the wide side up and drop the cartridge into the locking ring (with the small end down).

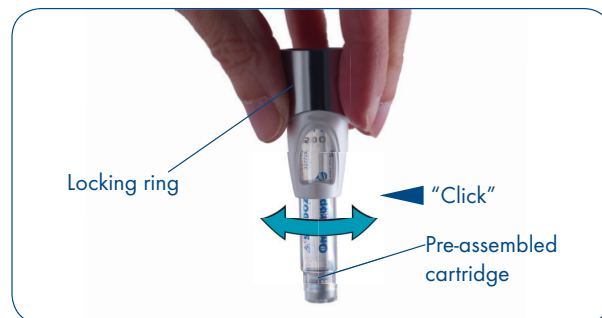

- Gently turn the cartridge until it drops down into the locking ring. Then pull the cartridge down until you hear it click into the locking ring. If the cartridge does not turn, it is already in the right position.

Attach locking ring to SurePal[®] body

- Screw the locking ring and cartridge back onto your SurePal[®]. Make sure the half circle on the locking ring and the dot on the SurePal[®] line up.

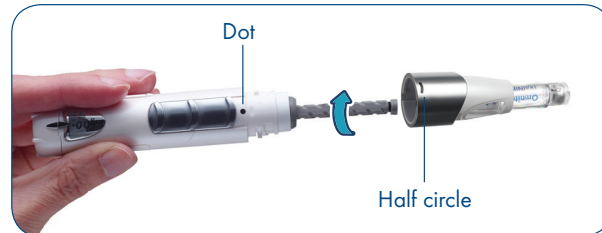

2. How to do a daily injection

Attach a needle


- Wipe the grey end of the cartridge with disinfectant (alcohol swab).
- Remove the needle paper seal from the needle and screw the needle onto your SurePal®. Take the outer cap off the needle, and keep it for later.



- Attach the needle hider (optional).

Turn the dose knob

- Turn the dose knob at the end of your SurePal® **gently**, until it stops. Your pre-set dose is now loaded. Check the dose window to make sure the dose is correct. In this example 2.4 mg has been set.
 - **Note:** if the dose in the window is not correct, unscrew the locking ring with cartridge. Then press down on the sliding button until “0.0” shows in the dose window. If necessary, then pre-set the dose as shown in section 1, and re-dial the dose knob.

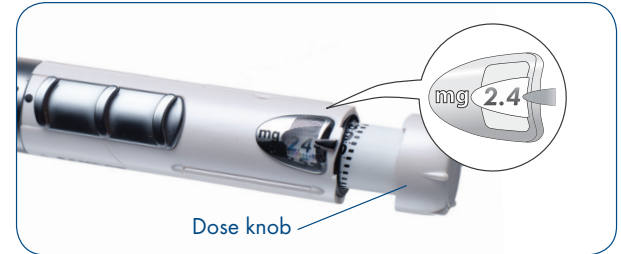

Insert the needle

- Select an injection site (e.g. thigh or belly), and clean it with disinfectant (alcohol swab). Rotate your injection site from time to time.
- Remove the inner needle cap and discard it. If you're using the needle hider, pull the sliding sleeve of the needle hider in a little bit so you can remove the inner cap. (Make sure your fingers are not covering the hole in the needle hider.)
- Pinch a fold of skin with one hand, hold your SurePal® upright, and insert the needle into the fold of skin. If you are using the needle hider, insert SurePal® into the fold of the skin, until the sliding sleeve is pushed all the way in. Release the skin when you push down on the sliding button.

Inject your dose

- Slide the sliding button down and hold. You will see the dose knob at the end of your SurePal® turning around.
- Hold the sliding button down for 10 seconds after the dose knob stops turning. Then remove SurePal®.
- Check that the dose window reads 0.0.

Note: if the dose does not show 0.0, change the empty cartridge and complete the injection with a new cartridge and new needle.
(See step 4: *Change empty cartridge during injection.*)



3. How to store your SurePal® Pen

Remove needle

- Attach the outer needle cap.
 - **Note:** if you used the needle hider, take the needle hider off before attaching the outer needle cap.
- Unscrew the covered needle, while holding the locking ring.
- Immediately discard the needle in a sharps box.
 - **Note:** needles must always be discarded in a special sharps box that you may get from your pharmacy or healthcare provider.



Store SurePal®

- Leave the cartridge in your SurePal®.
- Put the cap back on your SurePal®.
- Put your SurePal® in its hard case.
- Put the hard case in the refrigerator.

4. How to change an empty pre-assembled cartridge

Change empty cartridge after injection

- Remove and discard the needle in your sharps box.
- Unscrew the locking ring, and discard the empty cartridge.
- Put a new cartridge into the locking ring, and attach the locking ring to your SurePal®. (See step 1: *Attach the pre-assembled cartridge*).

Change empty cartridge during injection

If the dose knob stopped turning before the dose window reached 0.0, follow these instructions to complete your dose:

- Remove and discard the needle.
 - Unscrew the locking ring and discard the empty cartridge.
 - **Do not dial the dose knob or touch the sliding button!**
 - **Note:** your SurePal® knows how much is left in your dose, and the dose knob is already set to inject the right amount.
 - Put a new cartridge into the locking ring and attach the locking ring to your SurePal®. (See step 1: *Attach the pre-assembled cartridge*). Disinfect the grey rubber membrane of the cartridge with an alcohol swab before you attach a new needle. Then attach a needle, insert the needle and inject your dose (See step 2). Remember to hold the sliding button down for 10 seconds after the dose knob stops turning.
 - Store your SurePal® (see step 3).
-

OMNITROPE (somatropin) 5mg, 10mg and 15mg solution for injection administered using the SurePal pen is a fully funded (subject to Special Authority criteria, please refer to www.pharmac.govt.nz) **prescription medicine** used to promote the growth of long bones in children with reduced height due to lower than normal levels of growth hormone, to treat reduced growth in children with Turner syndrome and in children prior to the onset of puberty who have long-term kidney problems causing growth disturbance, and in children with Prader-Willi syndrome (PWS). In adults, it is used to treat those who have low levels of growth hormone and those who have PWS to improve their body composition.

Ask your doctor to see if OMNITROPE is right for you. You may need to pay for this medicine and normal doctor's fees apply. Use strictly as directed.

Omnitrope has risks and benefits.

Do not use OMNITROPE if you have an allergy to somatropin, the active ingredient or any of the other ingredients: monobasic sodium phosphate; dibasic sodium phosphate heptahydrate; poloxamer; mannitol; benzyl alcohol; phosphoric acid; sodium hydroxide; water for injection.

Do not use OMNITROPE if you are currently being treated for cancer, have active tumours or evidence of tumour growth; your bones are unable to grow anymore (closed epiphyses if you are a child); you have a serious injury or illness or surgical procedures requiring intensive therapy; you have Prader-Willi syndrome and are severely overweight or have difficulty breathing.

Do not give OMNITROPE to a child under the age of 3 years.

Tell your doctor if you have or have ever had diabetes, cancer, a kidney transplant, thyroid problems, are pregnant, planning to become pregnant or are breastfeeding, any disorder relating to hormones including adrenocortical insufficiency (also known as ACTH deficiency), are taking any other medicines.

Most common side-effects: a reaction near the injection site; limping; swelling of the arms or legs; painful joints; muscle aches, pains or stiffness; tingling; or numbness in your hands or feet. If symptoms continue or you experience other side effects or would like further information, please talk to your healthcare professional or visit www.medsafe.govt.nz for the Consumer Medicines Information.

OMNITROPE and SurePal are registered trademarks of Sandoz AG. Sandoz New Zealand Limited, Auckland.

SANDOZ

Sandoz New Zealand Limited, Auckland. Phone: 0800 726 369

MLR-0006-NZ February 2025 TAPS NP22283